

High quality



Safety

Reusable Laproscopic Instruments

Long-lasting performance. Numed offers an extensive selection of reusable detachable monopolar laparoscopic instruments available as well as in standard and bariatric lengths.



Ergonomic

Laprosopic Instruments Ending Types


Variety of jaw patterns of reusable inserts of scissors, graspers and dissectors makes the most complicated minimally invasive procedures possible. These include such specialities like General, Bariatric, Colorectal, Urological and Gynaecological surgery.



Durable



Excellent clinical performance and user experience

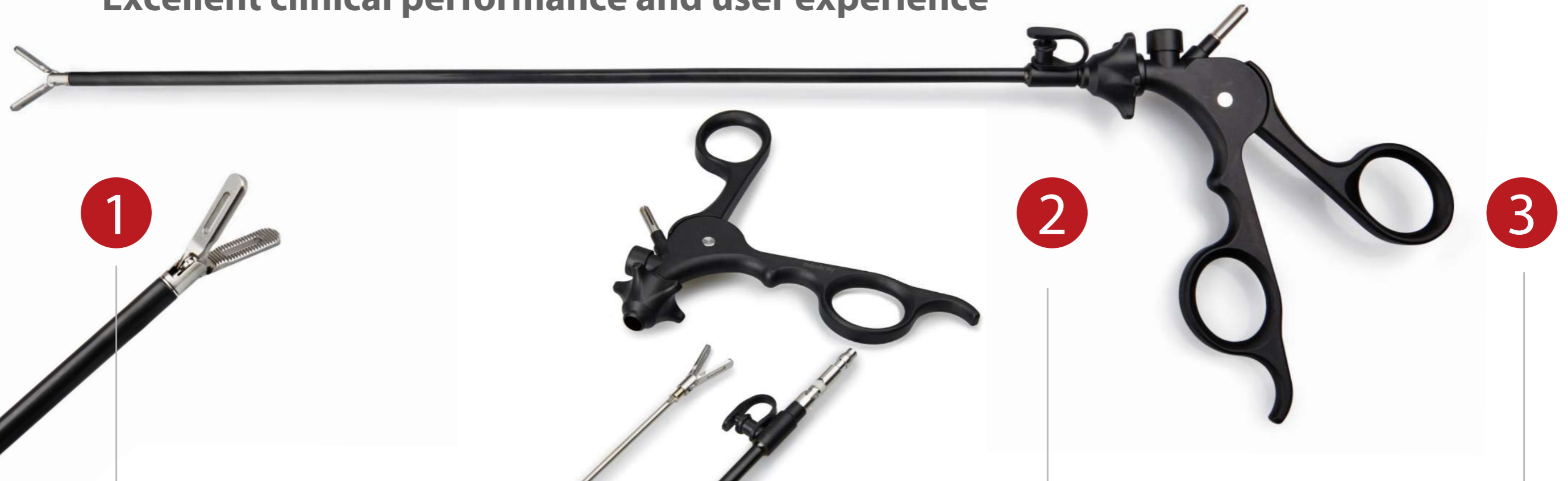

1

2

3

Jaws design



- Shorter metal conjunction effectively reduces risk of cautery damage ;
- Smooth hinged design reduces tissue and body fluid embedding.



- Refined jaws provides a strong grasping.

Modular design

- 3-part dismantling structure improves cleaning and sterilization efficiency and reduces the cost of replacement.

PEEK insulating sheath


- PEEK super insulation materials are wear-resistant and can withstand a voltage of 4000V.

Ergonomic design



- Thumb rest with extra-wide molding provides a more comfortable grip especially for prolonged surgeries.
- Light weight.

Simple assembly & disassembly

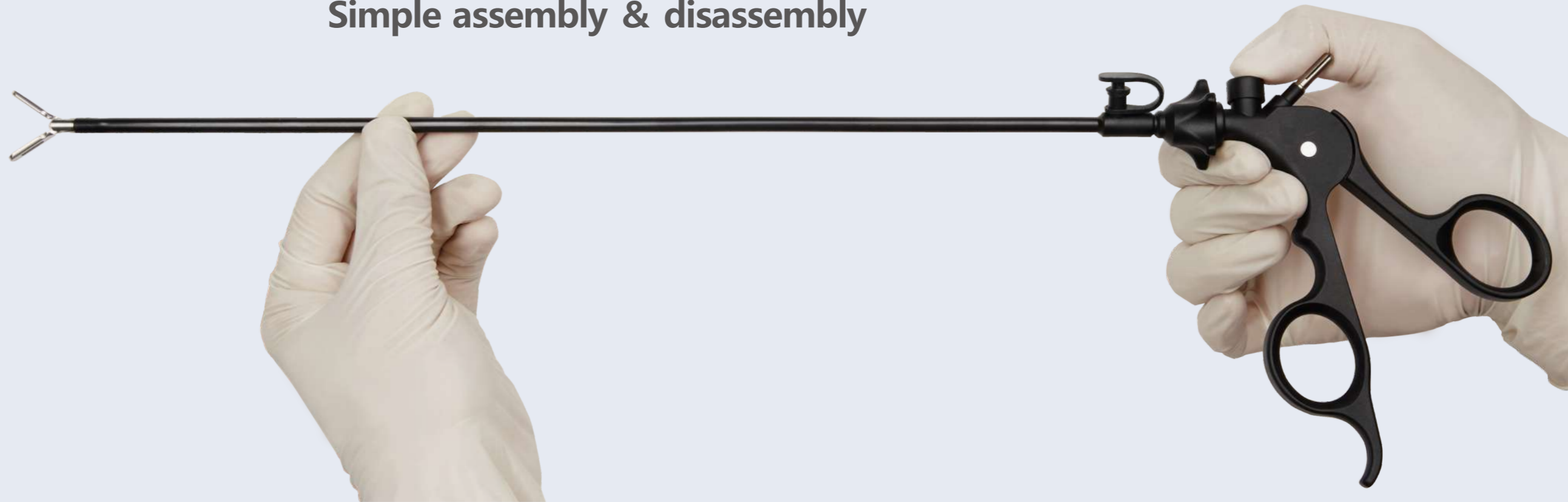

Assembly Process



- Put the insert in the shaft



- Rotate to lock the insert in the shaft



- Put the shaft into the handle



- Click and go !

Disassembly process



- Press the button



- Pull the shaft lever out of the handle



- Rotate the insert



- Pull the insert out

Other instruments

Endo Clip Applier




- Modular 3-part dismantling instruments for maximum hygiene;
- Handle with protection ratchet effectively prevents mis-clamping;
- All metal parts are crafted by CNC machines with high precision.

Bipolar forceps



- Special insulating coating effectively improves the durability;
- Perfect combination of reliable coagulation and strong grasping .

Needle Holder

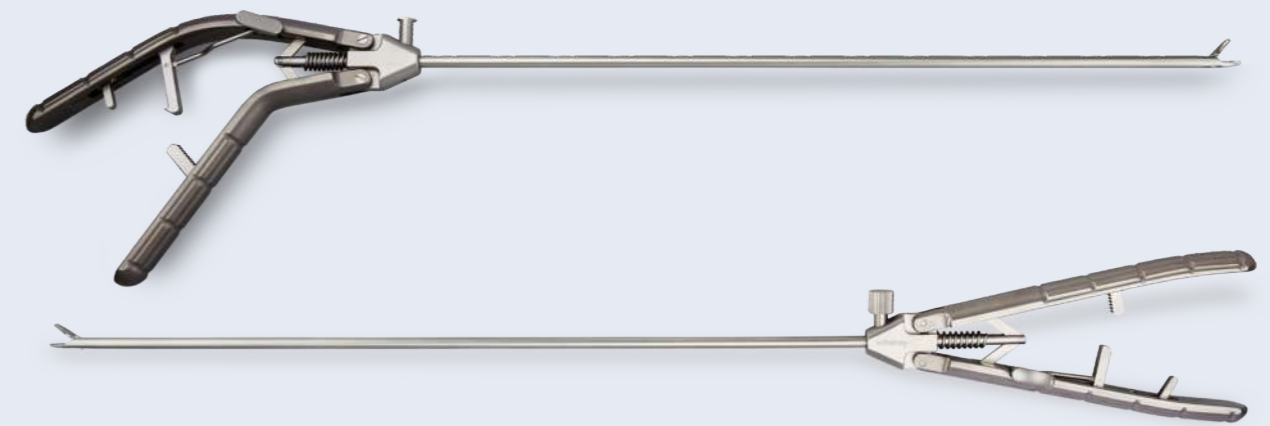

- Jaw welded with tungsten steel for stable holding and longer durability;
- Ergonomic handle and ratchet design for comfortable grasping and easy operation;
- Ultralight titanium alloy handle for lighter weight.

Scissors



One-step crafted blades



- Higher hardness;
- Smoother cutting;
- Longer life .

Monopolar Dissecting Forceps










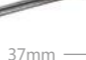

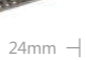

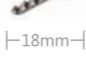

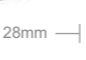

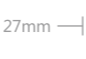


	
5mmX36cm	222-5Y301
10mmX36cm	222-0Y301
	

	Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)
	 Straight Dissecting and Grasping Forceps, Atraumatic	•		36cm	222-5Y3043	222-5S3043	222-5T3043
	 Kelly Dissecting and Grasping Forceps, Maryland	•		36cm	222-5Y3013	222-5S3013	222-5T3013
	 Kelly Dissecting and Grasping Forceps, Long, Maryland	•		36cm	222-5Y3023	222-5S3023	222-5T3023
	 Dissecting and Grasping Forceps, Atraumatic		•	36cm	222-0Y3103	222-0S3103	222-0T3103
	 Dissecting and Grasping Forceps, Right-angled	•		36cm	222-5Y3123	222-5S3123	222-5T3123
	 Dissecting and Grasping Forceps, Right-angled		•	36cm	222-0Y3133	222-0S3133	222-0T3133











Monopolar Grasping Forceps

	Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)
	 Manhes Grasping Forceps, "Tiger-jaws", 2x4 teeth	•		36cm	222-5Y3363	222-5S3363	222-5T3363







	
5mmX36cm	222-5Y301
10mmX36cm	222-0Y301
	

	Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)
	 Grasping Forceps, 2x4 teeth	•		36cm	222-5Y3393	222-5S3393	222-5T3393
	 Claw Forceps, 2x3 teeth, Short		•	36cm	222-0Y3403	222-0S3403	222-0T3403
	 Claw Forceps, 2x3 teeth		•	36cm	222-0Y3413	222-0S3413	222-0T3413
	 Sawalhe Tissue Grasping Forceps		•	36cm	222-0Y3423	222-0S3423	222-0T3423
	 Bowel Grasper, Fenestrated	•		36cm	222-5Y3303	222-5S3303	222-5T3303
	 Atraumatic Grasping Forceps	•		36cm	222-5Y3633	222-5S3633	222-5T3633
	 Dissecting and Grasping Forceps, "Alligator jaws"	•		36cm	222-5Y3773	222-5S3773	222-5T3773
	 Grasping Forceps, Atraumatic	•		36cm	222-5Y3333	222-5S3333	222-5T3333
	 Croce-olmi Grasping Forceps, Atraumatic, Curved	•		36cm	222-5Y3453	222-5S3453	222-5T3453
	 Manhes Dissecting and Grasping Forceps, "Duckbill"	•		36cm	222-5Y3283	222-5S3283	222-5T3283


Monopolar Grasping Forceps

						 5mmX36cm 222-5Y301 10mmX36cm 222-0Y301  	
Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)	
 Grasping Forceps, Atraumatic, Spoon-shaped <small>19mm</small>	•		36cm	222-5Y3223	222-5S3223	222-5T3223	
 Vancaillie Oviduct Forceps <small>13mm</small>	•		36cm	222-5Y3583	222-5S3583	222-5T3583	
 Vancaillie Adhesion Forceps, One Jaw Fenestrated <small>15mm</small>	•		36cm	222-5Y3493	222-5S3493	222-5T3493	
 Babcock Grasping Forceps <small>18mm</small>	•		36cm	222-5Y3553	222-5S3553	222-5T3553	
 Babcock Grasping Forceps <small>34mm</small>		•	36cm	222-0Y3573	222-0S3573	222-0T3573	
 Grasping Forceps, W Specially Fine Atraumatic Serration <small>26mm</small>	•		36cm	222-5Y3473	222-5S3473	222-5T3473	
 Bowel Grasper, Fenestrated, Short <small>18mm</small>	•		36cm	222-5Y3213	222-5S3213	222-5T3213	























Other Monopolar Forceps

						 5mmX36cm 222-5Y301 10mmX36cm 222-0Y301  	
Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)	
 Manhes Biopsy Forceps <small>15mm</small>	•		36cm	222-5Y3903	222-5S3903	222-5T3903	
 Spoon Forceps <small>32mm</small>		•	36cm	222-0Y4883	222-0S4883	222-0T4883	
 Debaquey Grasping Forceps, Curved and Slender Jaws <small>31mm</small>	•		36cm	222-5Y3813	222-5S3813	222-5T3813	















Monopolar Scissors

Product	5 mm	10 mm	Length	Insert Code	222-550 (S)	222-551 (T)
 Scissors, Spoon-shaped Blades, Serrated, Curved <small>20mm</small>	•		36cm	222-5Y3012	222-5S3012	


Needle Holder

	Product	5 mm	10 mm	Length	Code	
	V Handle Straight Needle Holder	•		33cm	221-56416	
	V Handle Left Curved Needle Holder	•		33cm	221-56426	
	Pistol Handle Straight Needle Holder	•		33cm	221-56516	
	Pistol Handle Left Curved Needle Holder	•		33cm	221-56526	
	V Handle Straight Needle Holder	•		33cm	221-56616	
	V Handle Left Curved Needle Holder	•		33cm	221-56626	
	V Handle Right Curved Needle Holder	•		33cm	221-56636	
	V Handle Reposition Needle Holder	•		33cm	221-56646	
	Pistol Handle Straight Needle Holder	•		33cm	221-56716	
	Pistol Handle Left Curved Needle Holder	•		33cm	221-56726	
	Pistol Handle Right Curved Needle Holder	•		33cm	221-56736	
	Pistol Handle Needle Holder with Reposition Function	•		33cm	221-56746	









Endo Clip Applier

	Product	5 mm	10 mm	Length	Code	
	Endo Applier for Medium Titanium Clip, Single Action		•	36cm	221-05176	
	Endo Applier for Medium Large Titanium Clip, Single Action		•	36cm	221-05186	
	Endo Applier for Large Titanium Clip Single Action		•	36cm	221-05196	
	Endo Applier for Small Titanium Clip, Single action	•		33cm	221-55116	
	Endo Applier for Medium Titanium Clip, Single action		•	33cm	221-05126	
	Endo Applier for Medium Large Titanium Clip, Single action		•	33cm	221-05136	
	Endoscopic Titanium Clip Applier for Large Titanium Clip single action Jaw		•	33cm	221-05146	















Open Surgery Applier

	Product	Length	Code
	Open Surgery Applier for Small Titanium Clip	18cm	221-01516
	Open Surgery Applier for Small Titanium Clip	28cm	221-00516
	Open Surgery Applier for Medium Titanium Clip	23cm	221-01526
	Open Surgery Applier for Medium Titanium Clip	28cm	221-00526
	Open Surgery Applier for Medium Large Titanium Clip	28cm	221-00536
	Open Surgery Applier for Medium Large Titanium Clip	28cm	221-00546

Trocar

	Product	Code	Description
	Magnetic Auto Shield Trocar	221-51215	Φ5x95mm
	Magnetic Auto Shield Trocar	221-01215	Φ10x95mm
	Pyramidal Magnetic Trocar	221-51115	Φ5x95mm
	Pyramidal Magnetic Trocar	221-01115	Φ10x95mm
	Auto Shield Trocar	221-51225	Φ5x95mm
	Auto Shield Trocar	221-01225	Φ10x95mm
	Auto Shield Trocar	221-81225	Φ12x95mm
	Magnetic Auto Shield Trocar with Screw	221-51235	Φ5x95mm
	Magnetic Auto Shield Trocar with Screw	221-01235	Φ10x95mm
	Pyramidal Magnetic Trocar with Screw	221-51135	Φ5x95mm
	Pyramidal Magnetic Trocar with Screw	221-01135	Φ10x95mm
	Auto shield Trocar with Screw	221-51245	Φ5x95mm
	Auto shield Trocar with Screw	221-01245	Φ10x95mm
	Auto shield Trocar with Screw	221-81245	Φ12x95mm
	Pyramidal Trocar with Screw	221-51145	Φ5x95mm
	Pyramidal Trocar with Screw	221-01145	Φ10x95mm
	Pyramidal Trocar with Screw	221-81145	Φ12x95mm
	Pyramidal Trocar	221-51125	Φ5x95mm
	Pyramidal Trocar	221-01125	Φ10x95mm
	Pyramidal Trocar	221-81125	Φ12x95mm

Trocar Reducer

	Product	Code	Description
	Trocar Reducer	221-61725	Φ15-Φ5x120mm
	Trocar Reducer	221-01715	Φ10-Φ5x20mm
	Trocar Reducer	221-01725	Φ10-Φ5x30mm
	Trocar Reducer	221-01735	Φ10-Φ5x50mm
	Trocar Reducer	221-01745	Φ10-Φ5x60mm
	Trocar Reducer	221-01755	Φ10-Φ5x95mm
	Trocar Reducer	221-01765	Φ10-Φ5x120mm
	Trocar Reducer	221-01775	Φ10-Φ5x160mm
	Mount-on Reducer	221-01815	Φ10-Φ5
	Cross Sealing Cap	221-01915	Φ10
	Cross Sealing Cap	221-51915	Φ5
	Sealing Cap	221-01925	Φ10
	Sealing Cap	221-51925	Φ5
	Sealing Cap	221-01935	Φ10-Φ5

Veress Needle



	Product	Code	Description
	Veress Needle	221-20127	Φ2.0x120mm
	Veress Needle	221-20147	Φ2.0x140mm
	Veress Needle	221-25127	Φ2.5x120mm
	Veress Needle	221-25147	Φ2.5x140mm

Monopolar Electrode





	Product	Code	Description
	Hook Tip Coagulation Electrode	222-51110	Φ5x330mm
	Ball Tip Coagulation Electrode	222-51220	Φ5x330mm
	Spatula Tip Coagulation Electrode	222-51440	Φ5x330mm
	Scalpal Tip Coagulation Electrode	222-52540	Φ5x330mm

Instruments with Metal Handle





Cholangiography Forceps



	Product	Code	Description	
	Cholangiography Forceps	221-55213	Φ5x330mm	

Bile Duct Knife & Dissecting Forceps



	Product	Code	Description	
	Bile Duct Knife	221-05239	Φ10x305mm	
	10mm Five Fans Dissecting Forceps	221-05913	Φ10x330mm	
	5mm Three Fans Dissecting Forceps	221-55913	Φ5x330mm	

Suction & Irrigation tube



Product	Code	Description
Suction & Irrigation Apparatus	221-55124	Ø5/Ø10x330mm

Monopolar Cable



Product	Code
Monopolar Cable	999-030