

Hwatime

Hwatime Medical. Make life more brilliant



To Establish a global leading brand



HT9

Plug-in Monitor



Beautiful Appearance Advanced Technology

- Plug-in modular design, various optional parameters, Direct display, easy to operate
- 17" colorful and clear TFT display, Touch screen optional
- Able to operate by keyboard or optional mouse. Supporting HD, VGA, BNC output
- Any integration for single parameter module or multi-parameter module
- 1 single modular slot and 1 double modular slot, Simultaneous display of at least 10 parameters
- 13 kinds of arrhythmia analysis, ECG waveforms of 7-lead displayed simultaneously
- 15 kinds of drug dose calculation

Clinical monitoring applications

Suitable for ICU/CCU, Operation Room, Post-operative Observation Ward, Emergency Room, & General Wards, etc.

Standard Parameters: ECG, RESP, NIBP, SPO2, PR, TEMP

Optional: IBP, ETCO2 modular, CO, EEG, AG, etc.




htwatime



Advanced Technology

- Anti-electrosurgical,
anti- defibrillation,
suitable for monitoring critical care.
- Reliable SpO2 technology,
good performance under the condition of motion and low perfusion.

Humanized Design, Perfect Expression

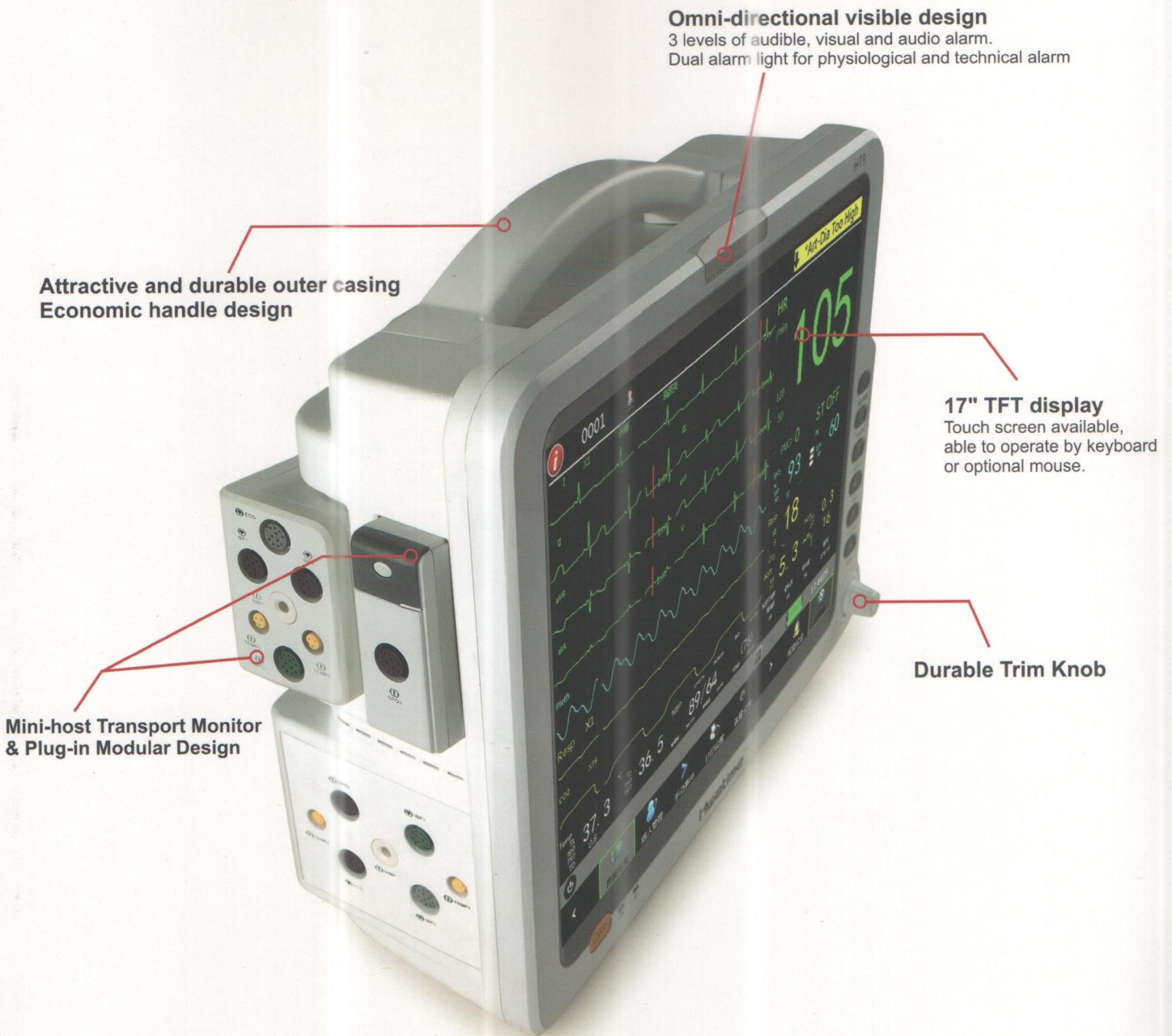

Hwatime

HT9 | Plug-in Monitor




Minihost Transport Monitor



Powerful software for complicated information processing

Supporting SD card, 8 USB interfaces

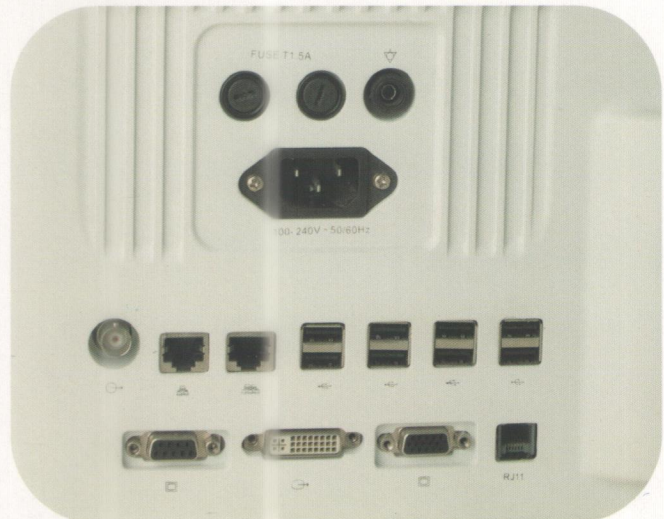

Hwatime



Hwatime



Multi interfaces display Comprehensive information presentation

Powerful software safe and stable


- Display of 1000 hours trend data and graphs
- Review of 48 hours full waveform
- Data store upon power failure
- Complicated information processing ability
- NIBP dual overpressure protection
- Wired/wireless connection
- Supporting SD card, 4 USB interfaces
- Various input
- Able to operate by keyboard or optional mouse



Standard Configuration

- ECG Cable: 1 piece
- TEMP Probe: 1 piece
- NIBP Cuff: 1 piece
- NIBP Extension Tube: 1 piece
- SPO2 Probe: 1 piece
- Power Cable: 1 piece
- ECG Electrodes: 20 pieces

HT9 | Plug-in Monitor

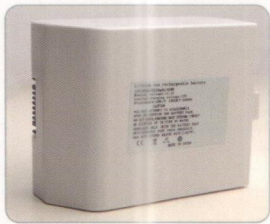

High Standard, High-quality Battery

Built-in detachable rechargeable lithium battery
Long battery life - 4.5 hours




Built-in Recorder

Able to print real-time waveforms


0 ST OFF
93 PR 60
120
50


8
InCO₂ 0.5 0.3
0.0
awRR 30 16
8

张压 平均压 时间
64 78 17:43:56

 17:44:09

 报警设定
 主菜单


SILENCE



PAUSE


FREEZE


NIBP



RECORD


MENU

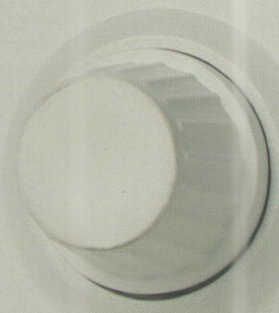

Durable Trim Knob

Omni-directional visible design,
3 levels of audible, visual, and audio alarm,
Dual alarm light for physiological and technical alarm.

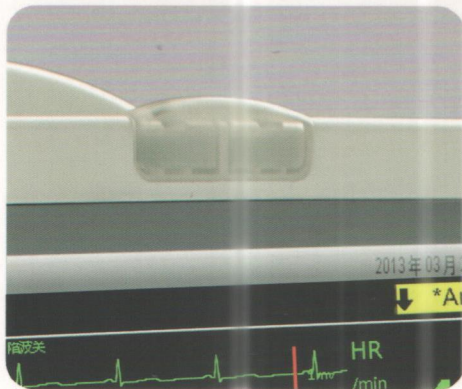


Plug-in Modular Design



17" TFT display

Touch screen optional, able to operate by keyboard or optional mouse

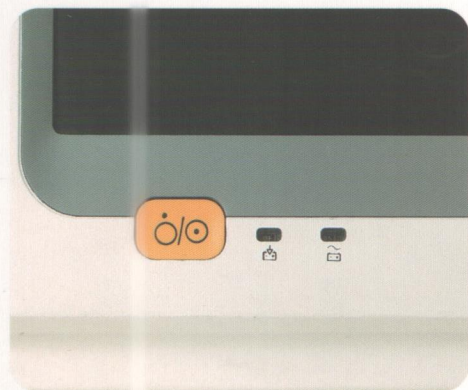




- Trends interface, easy to observe trends of each parameter
- Oxy CRG interface, more intuitive and centralized display of specific parameters



ECG Standard Interface

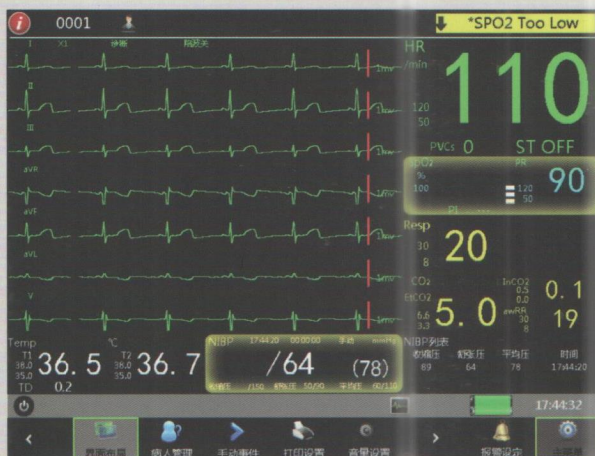

Full 7 lead ECG waveform Interface



HT Central Monitoring System
(Wired and wireless connections)



HT6



HT8



HT9

Hwatime

Hwatime Medical. Make life more brilliant

7/F,Block A,Bofook Jewelry Industry Park,No.33,
Bulan Rd,Longgang District,shenzhen city
Tel:0755-36678888 Fax:0755-85221277
E-mail:sales@hwatime.com Post Code:518112



Mobile 2-dimensional Code

For more information, please visit: www.hwatime.com Customer hotline: **4000 670 888**

- All rights are reserved by Hwatime. All technical specifications are subject to change without notice.
- See Operation Manual on prohibited contents or precautions!